

CARE CONCEPT
POZ' IN' FORM®
PHARMAQUEST INDUSTRIES

Poz'In'Form® cushions are designed to perform various support and comfort functions to support users in lying and at night time. The microfibers prevent from creep effects that are regularly observed in classical positioning and from early deterioration of the support. The various forms, sizes and models of Poz'In'Form® allow for a wide range of positioning possibilities.

- Poz'In'Form® elements and modules are entirely removable and washable. They are delivered with a Pharmatex or Lenzing FR® cover.
- The Pharmatex bi-elastic cover is made of polyurethane and is recommended for use in healthcare facilities. Pharmatex is soft to the touch and very comfortable next to the skin, while still easy to clean.
- The Lenzing FR® cover is made of viscose, coated with a waterproof polyurethane film. Lenzing FR® is a major innovation in positioning of the long-term bedridden, preventing maceration. Hygroscopic properties of viscose allow moisture absorption which is particularly important for people with thermoregulation issues. Moreover elasticity and flexibility of Lenzing FR® are particularly suitable for problems regularly observed in classical positioning.
- Lenzing spare covers can be purchased to cover Pharmatex pillows for a more absorbent feel.




• Pharmatex Cover (Wipe-down)



• Lenzing FR® Cover (Breathable, yet waterproof)

Poz'In'Form Delta

- This cushion is highly versatile, able to support users with a range of positioning needs. Under the knees, it provides support to the joint to maintain some flexion. Used in this way, the Delta also prevents users from sliding down in their beds which reduces friction and shear forces on the sacrum.
- It can be used in a semi-fowler position and helps to maintain positioning while relieving pressure in sensitive areas (sacrum, ischial tuberosities, heels) and shifting it to low-risk zones (calves, thighs and back).

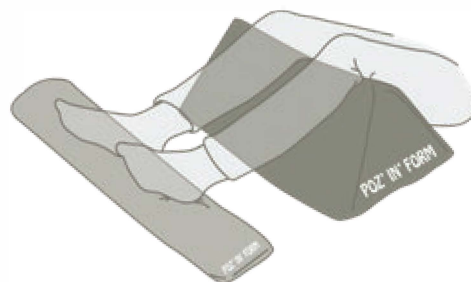

	With Pharmatex Cover	With Lenzing FR® Cover
Code	7009-0551-269	7009-0551-279
Length	56 cm	56 cm
Height	28 cm	28 cm
Depth	32 cm	32 cm

Poz'In'Form Universal

- Thanks to the three sizes, the universal cushion can be used in many different ways. The larger size can be useful for supporting knee flexion or moderate hip abduction when the user is in supine.
- The small cushions provide additional support to the user where needed, thus reducing the risk of abnormal postures or areas of high pressure.



	With Pharmatex Cover
Code	7009-1252-269
Length	55 cm
Depth	40 cm

	With Lenzing FR® Cover
Code	7009-1252-279
Length	55 cm
Depth	40 cm

Poz'In'Form Cylindrical

- Placed between the legs, the cylinder supports the upper hip and knee joints to be better aligned, reducing unwanted hip adduction and other unwanted changes to posture.
- In supine the cylinder can be applied behind the calves to alleviate pressure caused by contact between the users heels and the mattress surface.

	With Pharmatex Cover	With Lenzing FR® Cover
Code	7009-1501-269	7009-1501-279
Length	75 cm	75 cm
Diameter	18 cm	18 cm



Poz'In'Form Half-Moon

- Used in supine to provide support to the head, neck and upper limbs, including the shoulders.

	With Pharmatex Cover	With Lenzing FR® Cover
Code	7009-1851-269	7009-1851-279
Length	200 cm	200 cm
Width	30 cm	30 cm



Poz'In'Form Hand Wedge

- Cushion for the hands, wrists and forearms to provide greater surface area contact and distribute pressure more evenly.
- Can be particularly useful after certain surgical procedures where improved upper limb positioning is needed.

	With Lenzing FR® Cover
Code	7009-5411-279
Length	39 cm
Height	9 cm
Width	22 cm



Poz'In'Form Semi-Lateral Decubitus

- Used to provide body length support in semi side lying positions, from the top of the shoulders to the feet and ankles.

	With Pharmatex Cover	With Lenzing FR [®] Cover
Code	7009-3001-269	7009-3001-279
Length	185 cm	185 cm
Width	35 cm	35 cm



Poz'In'Form Ring & Half Ring

- In combination with the cylindrical cushion the half ring can provide medial thigh and knee support, preventing hip adduction and supporting a neutral position.
- Used for positioning the neck and head, providing a comprehensively supported position when used in combination with other cushions.

	With Pharmatex Cover	With Lenzing FR [®] Cover
Half Ring Code	7009-4051-269	7009-4051-279
Ring Code	7009-4081-269	7009-4081-279
Length	105 cm	135 cm
Width	22 cm	20 cm



Poz'In'Form Occiput

- Used to provide positioning and pressure relief to the occipital area in supine. Combination layers of microfibre and memory foam provide optimised pressure distribution.

	With Pharmatex Cover	With Lenzing FR [®] Cover
Code	7009-4751-269	7009-4751-279
Length	50 cm	50 cm
Height	9 cm	9 cm
Width	40 cm	40 cm



Poz'In'Form Heel Cushion

- Used at the bottom of the bed and fixed with an adjustable strap.
- Provide support behind the calves to raise the feet away from the mattress surface, thus reducing contact and associated areas of high pressure particularly at the heels.

	With Lenzing FR® Cover
Code	7009-5201-279
Length	85 cm
Diameter	25 cm



Poz'In'Form Heel Wedge

- Alternative cushion shape for supporting the calve and relieving pressure from the heel area.
- The large surface area of the cushion provides additional support for larger calve sizes or where more contact is needed to support the user.

	With Lenzing FR® Cover
Code	7009-5401-279
Length	48 cm
Height	12 cm
Width	30 cm

