
IC15 ACTIVEX™ MATTRESS






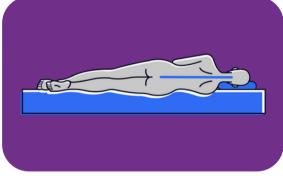

WITHOUT GOING TOO FAR

XP Support prevents you from bottoming out the mattress even on adjustable beds.

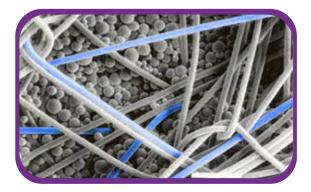

IMMERSE YOURSELF

ActiveX® provides optimal immersion and envelopment for long term comfort and pressure relief.




DONT SWEAT IT


Bluetec® and Vant® fabric technologies combine to keep your skin dry and cool.



BED CODES & SIZES

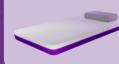

LONG SINGLE Code: IC15LS Mattress Size: 2030 x 900

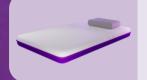

KING SINGLE Code: IC15KS Mattress Size: 2030 x 1070

LONG DOUBLE Code: IC15LD Mattress Size: 2030 x 1350

QUEEN Code: IC15Q Mattress Size: 2030 x 1520



SPLIT QUEEN Code: IC15SQ Mattress Size: 2030 x 750 (per side)

CARE & CLEANING

To launder the cover in a washing machine, set the main wash to a minimum of 65°C for no less than 10 minutes. The cover can be tumble dried on low heat and the temperature must not exceed 60°C. The foam can be autoclaved at 134°C.

SPECIFICATIONS

ActiveX[™] pressure relieving material Vant ventilated sides and base fabric Bluetec temperature regulating cover Option Zerotec waterproof cover (see pg. 22) Flip free design Zip off machine washable covers 45 - 300kg user weight Standard cover MVTR: NA. Medical Cover MVTR (g/24h·m2): 1635.0 Average Medical Cover BS7175: Crib 5 Waterlow ≤20

Feel Rating Mattress Heigh Risk Category User Weight Warranty Firm 150mm High Risk 45kg - 300kg 10 Year Warranty

PRESSURE CARE

Risk Category High

This rating is based on our average peak pressure of <24mmHg. For some high risk clients loare Zerotec covers are recommended

Icare risk category ratings are based on average peak pressure, inerapist recommendation, case studies, industry comparison and the Waterlow scale. For some high risk clients Icare Zerotec covers are recommended. The best way to assess a client on a pressure reducing surface is to monitor pressure points with visual checks over a period of 5 days.

backtosleep.com.au