


IC25 ACTIVEX™ MATTRESS



The IC25 is our most immersive mattress for users who prefer a soft feel. Suitable for light weight clients, those with bony prominences or pediatrics for pressure injury prevention and treatment.

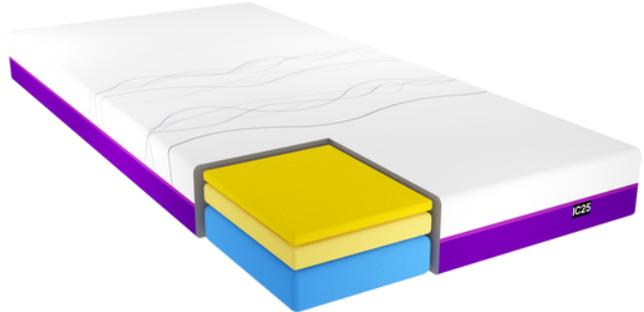

RISK CATEGORY:

HIGH RISK


COMFORT CATEGORY:

SOFT FEEL

KEY FEATURES



PRESSURE
REDUCING

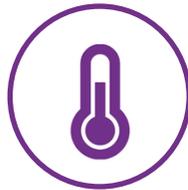

TEMPERATURE
REGULATING




BREATHABLE



ANTI SHEAR



MOISTURE
WICKING

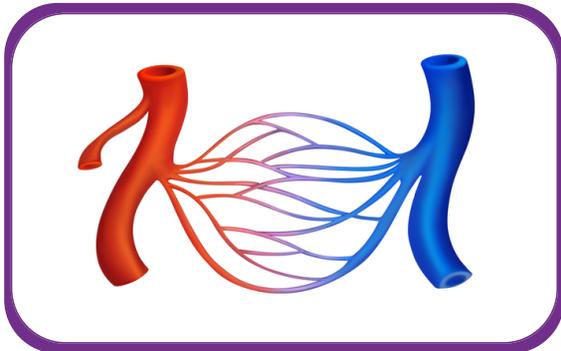


CIRCULATION

For thin/lightweight users, softer feel improves circulation.

FULL IMMERSION

Four layer construction maximises body contact reducing average peak pressure.

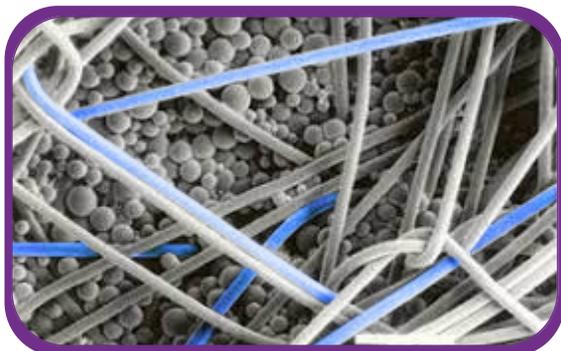

BLUETEC® AND VANT® FABRIC

Regulate surface temperature and help prevent peaks/spikes in temperature.


MORE FEATURES



10 Year
Warranty



200kg weight
rating

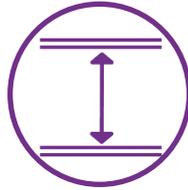

200mm mattress
height



Compatible with
Zerotec waterproof
covers



Flip free design



Zip off washable
covers



Multi layer
immersion



Anti Slip Base
Layer




Designed for use on
adjustable beds




20+ waterlow
rating


BED CODES & SIZES

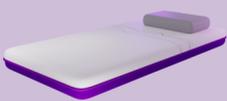

SINGLE
Code: ICXLS
Mattress Size: 1900 x 900

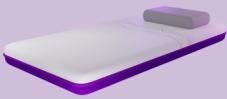


LONG SINGLE
Code: ICXLLS
Mattress Size: 2030 x 900

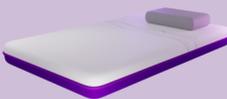


KING SINGLE
Code: ICXLKS
Mattress Size: 2030 x 1070




DOUBLE
Code: ICXLD
Mattress Size: 1900 x 1350

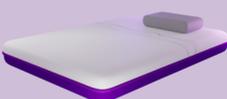

LONG DOUBLE
Code: ICXLLD
Mattress Size: 2030 x 1350

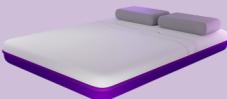

QUEEN
Code: ICXLQ
Mattress Size: 2030 x 1520

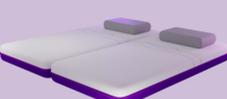

SPLIT QUEEN
Code: ICXLSQ
Mattress Size: 2030 x 750 (per side)

SPECIFICATIONS

ActiveX™ pressure relieving material
Vant ventilated sides and base fabric
Bluetec temperature regulating cover
Option Zerotec waterproof cover (see pg. 22)
Flip free design
Zip off machine washable covers
30kg - 200kg user weight
Standard cover MVTR: NA.
Medical Cover MVTR (g/24h·m2): 1635.0 Average
Medical Cover BS7175: Crib 5
Waterlow ≤20

Feel Rating	Soft
Mattress Height	225mm
Risk Category	High Risk
User Weight	30kg - 200kg
Warranty	10 Year Warranty

PRESSURE CARE

Risk Category High

This rating is based on our average peak pressure of <22mmHg. For some high risk clients Icare Zerotec covers are recommended.

Icare risk category ratings are based on average peak pressure, Therapist recommendation, case studies, industry comparison and the Waterlow scale. For some high risk clients Icare Zerotec covers are recommended. The best way to assess a client on a pressure reducing surface is to monitor pressure points with visual checks over a period of 5 days.

CARE & CLEANING

To launder the cover in a washing machine, set the main wash to a minimum of 65°C for no less than 10 minutes. The cover can be tumble dried on low heat and the temperature must not exceed 60°C. The foam can be autoclaved at 134°C.

For Sales: 1300 848 077

backtosleep.com.au